

QUICK, EASY AND ECONOMICAL!

PORTER Scavenger Retrofit Kits

UPGRADE your present scavenger system to PORTER's high level of efficiency and safety

All kits include the Porter Scavenger with In-Line Vacuum (specify "A" for Adult Hood or "P" for Peds Hood), Three Liter Breathing Bag, Flowmeter Adapters and Bag Tee where necessary.

**PORTER
UNIVERSAL
RETROFIT KIT
#5800-1**

PORTER AVS 5000 DOES NOT RETROFIT TO OTHER MANUFACTURERS' FLOWMETERS

THIS ONE KIT RETROFITS ALL THESE FLOWMETERS!

- **MATRIX**
Quantiflex, MDM
 - **OHMEDA**
 - **FRASER HARLAKE**
 - **FRASER SWEATMAN**
 - **MDT-MCKESSON**
Analog Four/Compalor
 - **PARKELL-SEDATRON**
(not pictured)
 - **BELMED**
(not pictured)
 - **ACCUTRON**
(not pictured)
- Also replaces: Porter Clean-Air™ Rubber Goods Non-Scavenging System (Mask with plastic sweeps)

| Porter Part # | Retrofits |
|---------------|---|
| # 5800-2 | Coastal NCG Chemetron |
| # 5800-4 | Hampton - Simplaire Micro |
| # 5800-5 | McKesson - Analog I & II Starflite Denetron |
| # 5800-6 | McKesson - Analog 3C |
| # 5800-9 | Porter Clean-Air™ Rubber Goods & Canister Scavenging System |

NOTE: It is not recommended to retrofit Porter Scavenger Rubber Goods to Adec, Veriflo or Spica machines. Contact your authorized Porter Dealer for Trade-in offer.

Upgrade & Service Your Flowmeter

Flowmeters should be calibrated and have maintenance service every 2 years.

There are literally thousands of flowmeters in dental operatories that have not been checked in over a decade. Many of the manufacturers of many of these units are now out of business.

Porter Instrument Company offers fast turn-around service on the following flowmeters: Porter Standard Serial No. 14000-25000; Porter MXR Serial No. 50000-80000.



MATRIX
Quantiflex, MDM



OHMEDA



FRASER HARLAKE



FRASER SWEATMAN



MDT-McKesson
Analog Four



BELMED



COASTAL



NCG CHEMETRON



HAMPTON
Simplaire Micro



McKesson
Analog I & II



STARFLITE
DENETRON



McKesson
Analog 3C

PORTER

NITROUS OXIDE SEDATION SYSTEMS



PORTER INSTRUMENT COMPANY, INC.



245 TOWNSHIP LINE RD. O.P.#0907 HFIELD PA440-0907 USA 15772-4000 / FAX 15723-1299

FM-791 Rev. E 08/03



New AVS Technology... Exclusively From Porter

Our unique patent approved design automatically activates the scavenger vacuum as soon as gas begins to flow from the flowmeter. As a result, it prevents Nitrous Oxide/Oxygen from being administered without scavenging. In accordance with NIOSH guidelines, the vacuum should always be "ON" while administering Nitrous Oxide. The AVS 5000 retrofits to all flowmeters manufactured by Porter Instrument. The complete system now includes the flowmeter, the AVS 5000, and the scavenger breathing circuit.

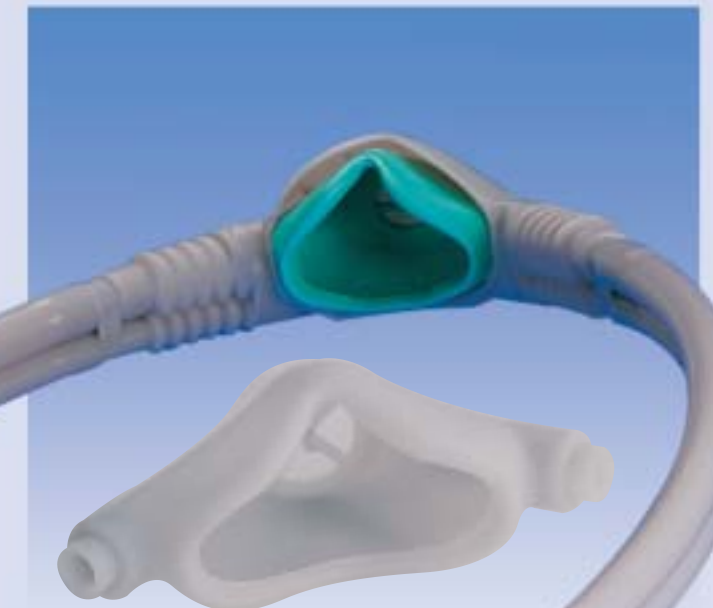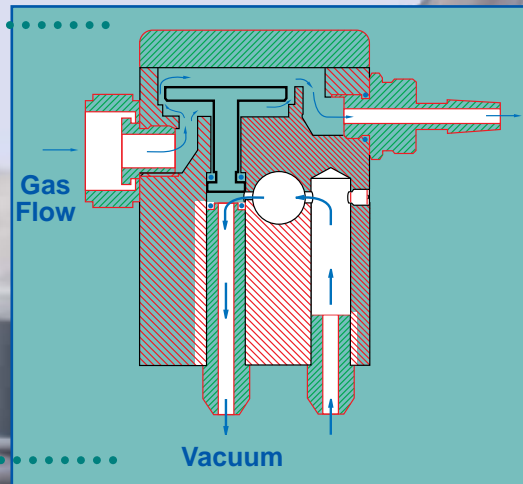

Gray Autoclavable



Green Disposable*



White Citrus Disposable*

When Gas Flows – Vacuum is "ON"... It's Automatic

When the procedure is finished, the vacuum is automatically shut down by the AVS. This is accomplished by turning off the flowmeter's Flow Control Knob.

Convenient Operation –

The vacuum level sight glass lets the user monitor the levels of vacuum and make adjustments to the flow rate when necessary. Porter also designed the system to reduce tubing clutter.

Worry-Free Scavenging System

- A flexible inner liner creates a tight seal around the nose for improved scavenging. Porter's unique double mask sets the industry standard for comfort and scavenging efficiency.
- Masks are autoclavable and available in Adult and Peds sizes. Disposable mask liners are also available for offices looking for added convenience.

Autoclavable or Disposable Liners –
***Disposable Liners are available in Unscented or Citrus scented.**